



EMBRACE NEST

Infant warmer



INTENSIVE CARES



Description

- Non-electric infant warmer
- For treating and preventing infant hypothermia in prehospital transport, intrahospital transport, and as an additional source of warming when incubators or radiant warmers are not available.
- Consists of a blanket, which has been specially designed to accommodate an infant, and a warming pack which is installed in the back zipper pocket of the blanket. Together, they provide instant warmth to assist in maintaining an infant's body temperature of 37°C.
- Generates warmth for at least 4 hours at ambient temperature greater than 25°C.

Measurements

Size: Medium

Embrace Nest dimensions: 52x25 cm 0,3 kg

Heater dimensions: 44x29x7 cm 2,5 kg

Warm pak dimensions: 38x22 cm 1,3 kg



WARMPAK
Maintains temperature of 37°C for at least 4 hours*
*duration of at least 4 hours at ambient temperature greater than 25°C
380x220 mm - 1.3Kg



Special Phase Change Material in WarmPak™ maintains a temperature of 37°C for at least 4 hours*



Portable and simple to use




Only 15 Watts of power consumption in ready mode



Does not require a constant supply of electricity



BABYWRAP
Reusable and easy to csanitize using disinfectants or soap wash
520x250 mm - 0.3Kg

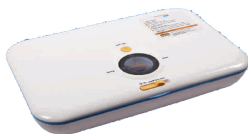


ACCUTEMP HEATER
Precision heating mechanism prevent overheating of WARMPAK
440x290x70 mm - 2.5Kg
230V, 50Hz



Easy to sanitize and reuse

Simple to use

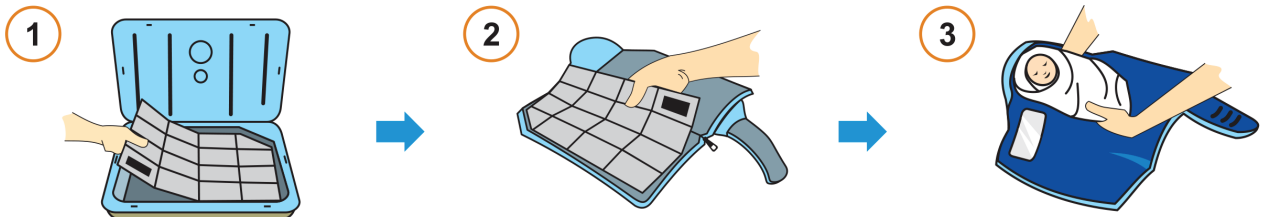

ESSE3 srl, Via Garibaldi 30
14022 Castelnuovo D.B. (AT)
Tel +39 011 99 27 706
Fax +39 011 99 27 506
e-mail esse3@chierinet.it
web: www.esse3.dreamgest.com

