



# ES-QCT Series

## Pathology Workstation



LABORATORY



## Description

Pathological sampling bench is widely used in hospital pathology department, pathology laboratory, etc. Reasonable ventilation system protects operator away from the harmful gas produced by the formalin during sampling. The hot & cold water system ensures it can adapt work in different climates.

### ADVANTAGE

1. Odor-remove function.
2. Single water sink, two-person operations.
3. Exhaust air system: High-power suction device, low noise.
4. Optional water heater is available, realizing hot & cold water adjustment.
5. The base cabinet is inward-concave type, convenient for sitting operation.
6. Optional pulverizer can completely smash the waste during test, to avoid blocking drainage pipes.
7. Automatic flushing device, slightly right-side inclined worktable, easy for waste discharge and cleaning.



### Configuration

UV Light; Fluorescent Light; Adjustable Spotlight; Stainless Steel Hanger; "UV, FAN, LIGHT, PULVERIZER" Control.




### Water Tap & Sink

Stainless steel; Whole-structure.



### Suction Inlet

Suction inlet adopts louver-type. The upper mouth structure ensures above blower suction, avoiding air back-streaming.

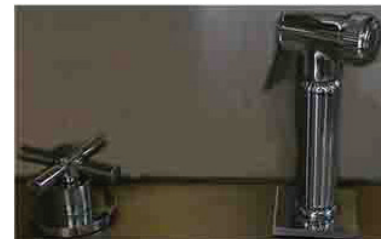

### Auto Flushing Design

Automatic flushing device, slightly right-side slope worktable, easy to waste discharge and cleaning.



Flex Magnifier Light(Optional)



### Pulverizer (Optional)

Completely smash the waste during test, to avoid blocking drainage pipes.



## Features

Model	ES-QCT-1000	ES-QCT-1500	ES-QCT-1800
External Size (W*D*H)	1000*750*1950mm	1500*750*1950mm	1800*750*1950mm
Working Area Size (W*D*H)	895*430*900mm	1395*430*900mm	1695*430*900mm
Work Surface Height	900mm		
Consumption	250W	300W	300W
Optional water heater	1500W		
Power Supply	AC220V±10%, 50/60Hz; 110V±10%, 60Hz		
Exhaust Airflow Volume	500-600m <sup>3</sup> /h		
Material	Main Body	304 stainless steel	
	Work Table	304 stainless steel	
	Standard Accessories	Spotlight; UV lamp; Fluorescent lamp; Exhaust system; 4m exhaust duct; Pipe strap; Auto flushing device; Stainless steel hanger; Adjustable hot and cold water tap; 180mm deep stainless steel water sink, funnel with filter. Active carbon filter	
	Optional Accessories	Pulverizer, Electrical water heater, Automatic induction hand dryer, Eye washer, Camera. Flex magnifier light, laptop bracket	
Package Size(W*D*H)mm	1160*890*2170	1660*890*2170	1960*890*2170
Gross Weight	250kg	300kg	340kg



## Accessory

### Standard Accessories

Spotlight	1pc	
UV lamp	1pc	
Fluorescent lamp		
	1pc	
Exhaust system	1pc	
Auto flushing device	1pc	
Stainless steel hanger	1pc	
Adjustable hot and cold water tap	1pc	
180mm deep stainless steel water sink, funnel with filter		1pc
Active carbon filter	1pc	

### Optional Accessory

Pulverizer  
 Electrical water heater  
 Automatic induction hand dryer  
 Eye washer  
 Camera



ESSE3 srl, Via Garibaldi 30  
 14022 Castelnuovo D.B. (AT)  
 Tel +39 011 99 27 706  
 Fax +39 011 99 27 506  
 e-mail [esse3@chierinet.it](mailto:esse3@chierinet.it)  
 web: [www.esse3.dreamgest.com](http://www.esse3.dreamgest.com)

