



PRESSURE REGULATORS FOR OXYGENTHERAPY

Mod. **FM**

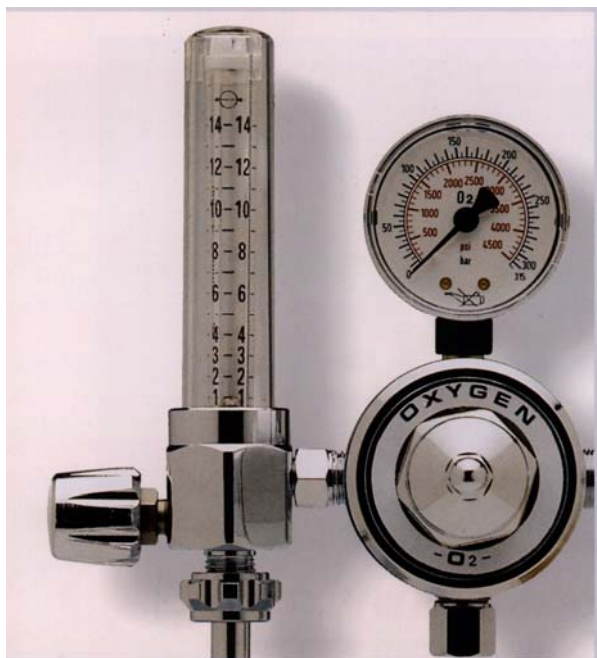

Other standard features, besides the ones already mentioned are

- Flowmeter type RS, 15 l/m' for O₂ at 3,5 bar, compensated pressure with outlet connection threaded ¼ ISO 3253 M;
- Humidifier type TR200 with inlet connection threaded ¼ ISO 3253 F. and outlet tubing connector Ø 6 ÷ 9.

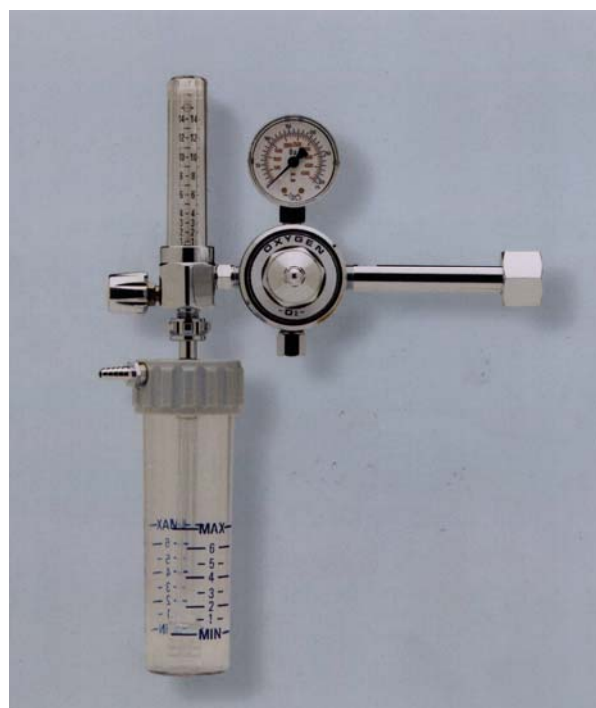

The pressure regulators for oxygen therapy type FM are suitable for the direct mounting on the cylinders, they are manufactured with the inlet connection for oxygen according to the different European standards and they are supplied complete with a compensated pressure flowmeter and a bubbling humidifier. Available with present outlet pressure calibration to check the cylinder contents, the oxygen therapy pressure regulators type FM have the body and the nipples made of chrome-plated brass, the plastic parts made of polycarbonate and autoclavable at 134°C, they are fitted with a pre-calibrated safety valve and they have a tubing connector on the outlet to patient.:



Laboratorio e Assistenza:

Via Garibaldi 30 - 14022 castelnuov DB (AT)
Tel: 039 - 011 99 27 706 Fax: 039 - 011 99 27 506
www.4ci.com E-mail info@4ci.com

