

Surgical INOX washbasin **ESMP-DUO®** Electronic control

Aseptic surgical sink with two washing stations
ESMP-DUO® allows washing up to the elbows
thanks to its wide and deep tank.

The enveloping shapes provide anti-splash
protection from the environment, preventing water
stagnation and the risk of bacterial development.

Without overflow, with protective wall backsplash
and rounded shapes with internal slope. Stainless
steel frame and HDPE supply hoses resistant to
chlorination.

Material: Satin polished 304 stainless steel

Departments concerned: Operating theater,
intensive care, nephrology, hemodialysis,
emergency box, surgeon preparation room.

Dimensions

L: 1400 mm
H: 1020 mm
P: 475 mm
Tank depth: 310 mm
Weight: 34.4 Kg




Control type

Electronic

Standard Equipment

- 2 IR cells triggered by presence detection
 - + Power supply of your choice by batteries or mains. Supplied 2 in 1
 - + System equipped with automatic purge (6h/12h/Off)
 - + Simplified adjustment on the front via remote control
 - + EV box easy to unclip for maintenance
- Removable straight spouts in chrome-plated brass
- 1 Ø40mm drain and siphon
- 1 Flow limiting valve
- Connection hoses
- Stainless steel wall mounting kit

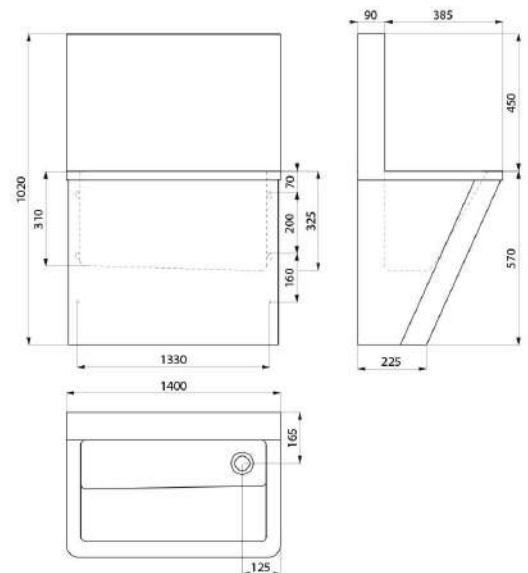

Optional accessories

- Thermostatic mixer + Anti-scald safety
- Mini pre-filtration kit with 1 or 5µ cartridge
- 40W ultraviolet kit
- Dry siphon without water retention
- Disinfectable hospital siphon
- Autoclavable detachable spout
- 0.2µ terminal filter for all germs, 14 or 31 days
- 60-day anti-legionella terminal filter
- 0.2µ terminal filter for all germs, 2-3-4 months
- Aluminium accessory support rail
- Manual or elbow mechanical soap/SHA dispenser
- Electronic soap/SHA dispenser
- Surgical nail brush dispenser
- Basket for hygiene accessories
- Paper dispenser